


Protevo GTE+ MCM

Convertible alternating pressure with microclimate management





Protevo GTE+ Microclimate is an advanced pressure redistribution mattress that provides flexible therapy to help you meet the wound care needs of your patient. Used without a control unit, Protevo GTE+ provides pressure redistribution through an open, self adjusting air cylinder system. Easily convert the surface with an optional control unit to provide microclimate management and alternating pressure.

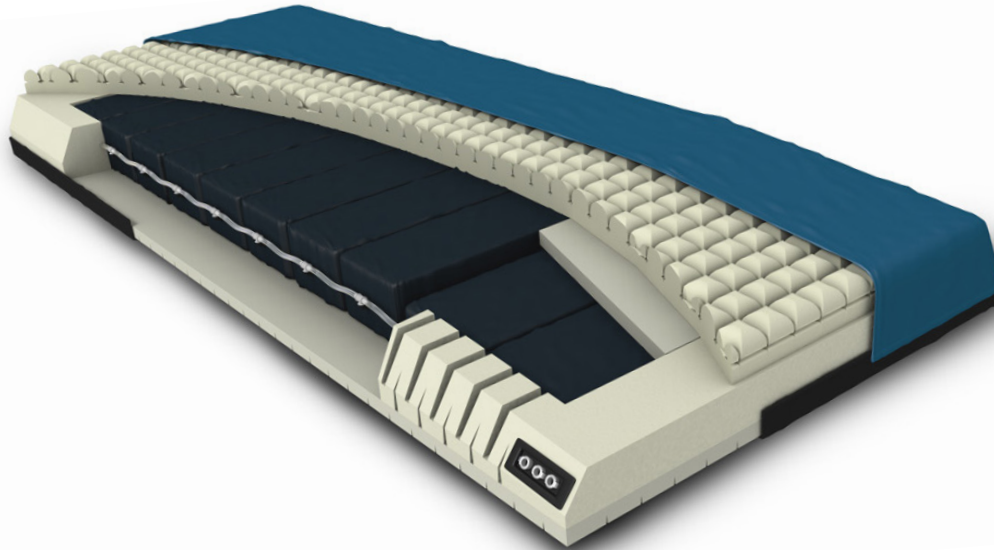

protevoTM GTE+ microclimate

- Advanced pressure redistribution mattress can be used with or without power
- Indicated for prevention and treatment of pressure injuries through Stage IV
- Open air system with self-adjusting air cells
- 4-way stretch cover and 3-dimensional sculpted envelopment layer provide optimal patient comfort
- 10 self adjusting air cells provide optimal pressure redistribution with dedicated air cells for the heels
- Accommodates patients weighing up to 550 lbs.
- Optional control unit provides microclimate management and alternating pressure with quiet operation - no noisy, constant adjustments of an all air mattress system

AT A GLANCE

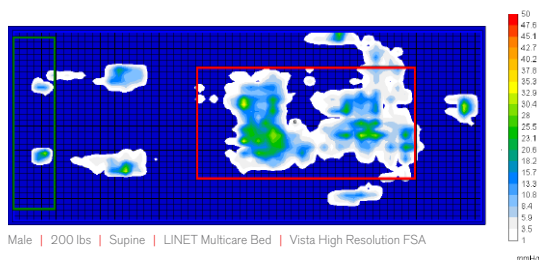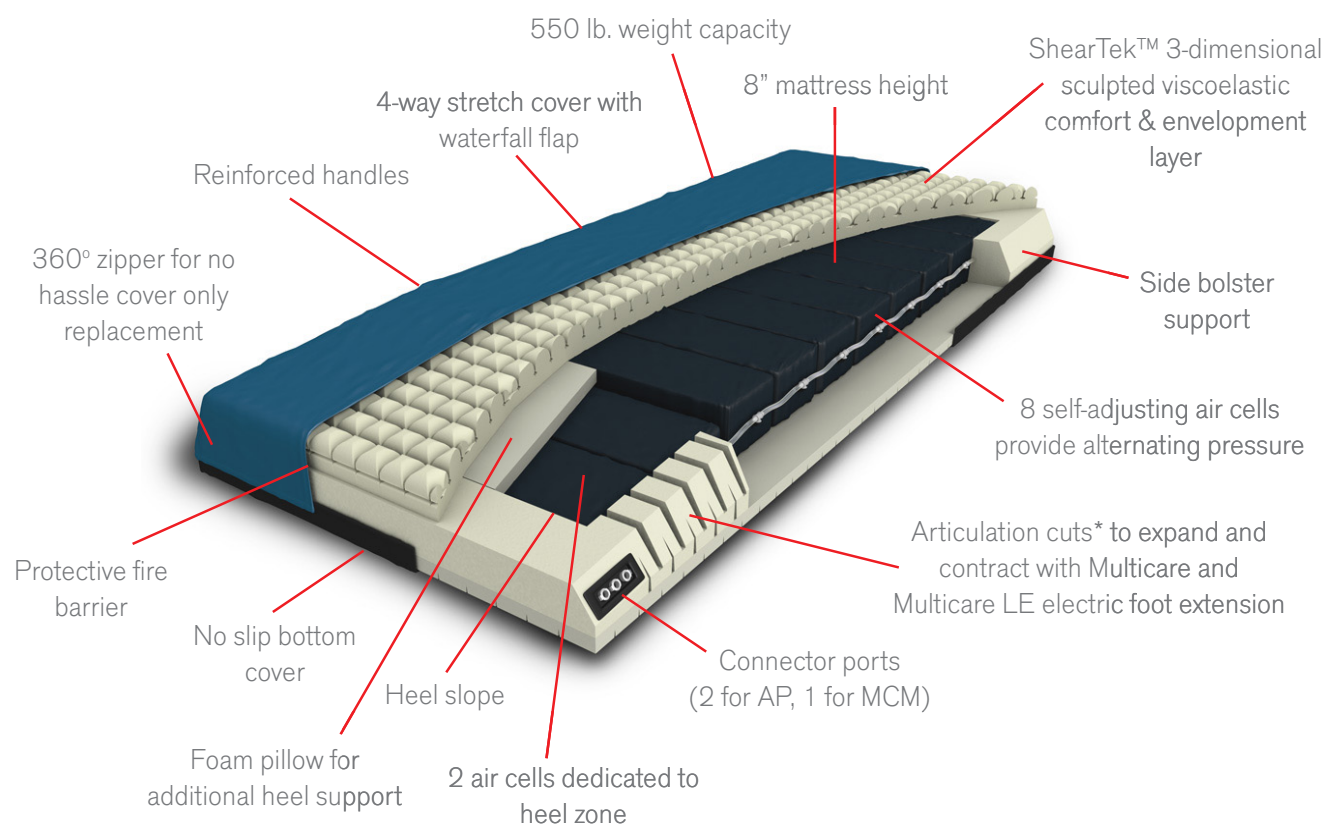

Sensing Area (in ²)	505.46	11.23	365.05
Coefficient of Variation (%)	83.70	76.08	80.63
Average Pressure (mmHg)	8.79	11.51	9.24
Maximum Pressure (mmHg)	36.48	32.17	36.48
Regional Distribution (%)	100.00	2.91	75.89

TECHNICAL SPECIFICATIONS

Mattress

L x W x D (flat deck)	84 x 35 x 8 in.
L x W x D (Multicare)	85.5 x 35 x 8 in.
Weight	26.5 lbs.
Therapeutic Weight Limit	550 lbs.
Flammability	16 CFR 1633
MVTR	500 g/m ² /24h
Mattress Warranty	5 yrs, prorated
Cover Warranty	2 yrs, prorated
Flammability**	Cal 117, Cal TB129, CPSC 16CFR1633, Boston 1X-11

Optional Control Unit

L x W x D	10.25 x 4.75 x 8.5 in.
Weight	5.5 lbs.
Airflow	10 lt./min.
AP Cycle Time	12 min.
Noise Level	35dBA
Warranty	2 yrs, non-prorated
Power Supply	120-127 VAC, 60 Hz
Operating Temperature Range	32° - 104° F
Relative Humidity Range	30-75%

* Articulation cuts only available on the Multicare and Multicare LE bed versions. **Protevo GTE+ meets or exceeds noted fire code standards.