

# ergozone<sup>TM</sup> iGEL

Gel Infused Advanced Pressure Redistribution Mattress



**Wound prevention. Superior Comfort.**

Ergozone<sup>TM</sup> iGel advanced pressure redistribution mattress is designed specifically for patient comfort and pressure ulcer prevention. Specialized gel technology, infused into a honeycomb-molded mattress core, provides a micro-zoning effect that supports and envelops the patient on an air suspension-like system. Surface Modification Technology works in combination with the infused gel to minimize pressure, maximize envelopment and increase patient satisfaction.


# Ergozone iGel – advanced pressure redistribution mattress



## Highlights

- Micro-zone pressure redistribution
- Top cover helps manage shear and friction
- Maintains ambient temperature
- Increased patient comfort and satisfaction
- Weighs over 50% less than other gel-based mattresses

The Ergozone iGel mattress can be used with any flat deck frame. When used with a LINET bed, the Ergozone iGel works in combination with the Ergoframe® design of the bed to help maintain patient position, neutralize shear and minimize pressure in the sacral area. This leads to improved patient comfort and easier and safer nursing care.

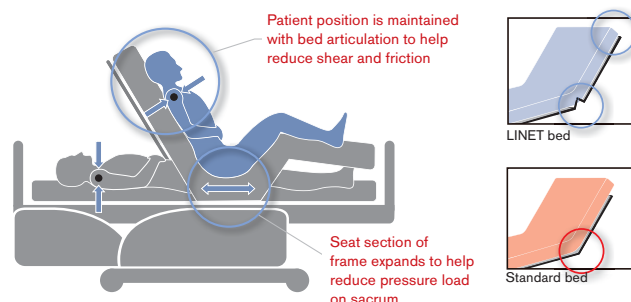

## SURFACE MODIFICATION

**TECHNOLOGY (SMT)** creates alterations in the mattress (by pattern, size, depth, spacing and location) to maximize pressure redistribution and protect the most vulnerable areas of the body.

**OPEN CELL TECHNOLOGY** removes body heat at a greater rate than conventional foams providing a cooler, drier experience for the patient.

- Greater air flow for temperature and moisture management
- Increased support retention over the life of the mattress

**INFUSED GEL TECHNOLOGY** works in combination with SMT to allow envelopment and support on a micro level.

- Enhanced patient comfort and satisfaction
- Maintains ambient temperature and provides a localized cooling effect on the skin

ergozone<sup>TM</sup> iGEL

**5-year Warranty**



**Environmentally Friendly Manufacturing Process** with natural ingredients for sustainability



**Enhanced Pressure Redistribution Zones** engineered for high risk areas including the head, shoulders, hips, legs and heels



**Envelopment Zones** strategically combine with body contouring zones to enhance pressure redistribution



**Patient Weight Support up to 500 lbs**



**Open Cell Visco-Elastic Foam** to allow greater air flow for temperature and moisture control



**4-Way Stretch Cover** moves in every direction to minimize shear and friction



**Gel Infused Surface** offers enhanced heat dissipation, moisture management, pressure redistribution and patient comfort

