

# AVE 2



# The journey begins

The innovative design of the AVE 2 brings a new level of comfort and efficiency to every phase of childbirth.



Attractive and functional design.

WUTAAATZOG



Easy and quick to clean.



AVE 2 is an ideal platform for safely performing labor exercises with a low height of 23.6 inches.

Flexible positioning options allows the bed to conform to the needs of the patient.




The maximum height of 41.3 inches decreases strain on caregivers and improves access for procedures.





# Benefits of AVE 2



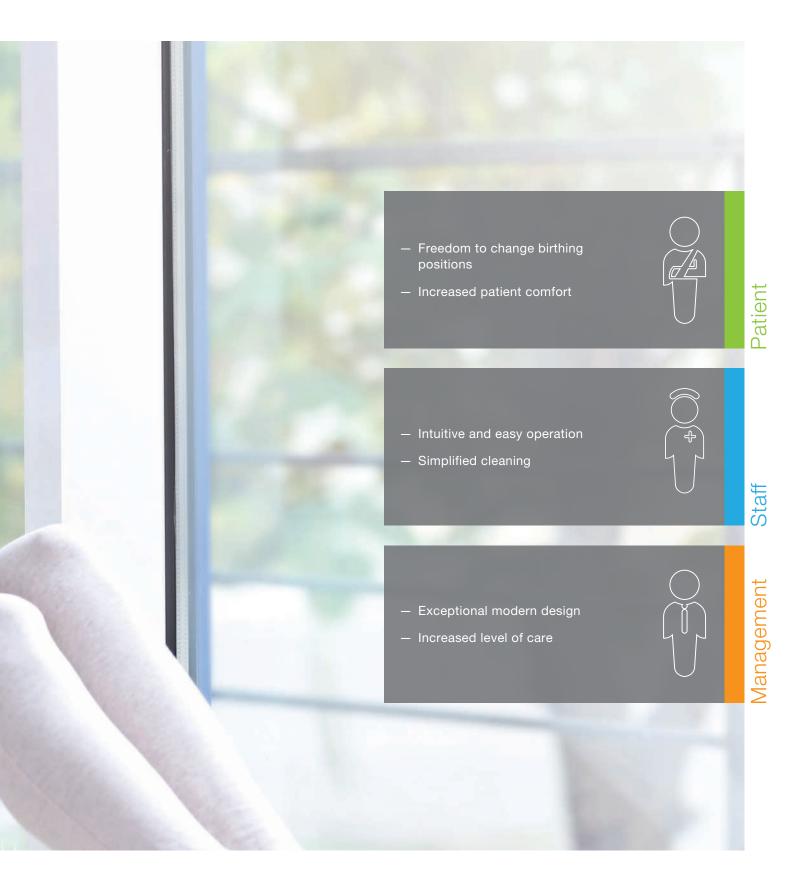

## Innovation

Simple swivel mechanism allows the foot section to retract and rotate away from the patient quickly and easily.

- Improve caregiver safety and help reduce injury by eliminating heavy foot section lifting.
- Reduce contamination and trip hazard from foot section never touching the floor.
- Facilitate rapid and efficient bed configuration.
- Enhance access to patient.



The innovative foot section allows for quick configuration and new positioning options for caregivers and patients.



When tilted, the section can be also used as a foot support.



# Simple operation

Rapid configuration with simple and intuitive controls improves efficiency and reduces caregiver frustration.



Quick adjustment of the leg rests into the required position.



Height adjustment for patients of all sizes.



Swing away for improved access.

### U-SHAPE

The unique oval shape improves perineal access for procedures and minimizes gap around fluid basin.





HAND GRIPS

Stable and comfortable grips.

## COLUMN CONCEPT

Innovative column design allows for both a low height for patient egress and greater height for improved caregiver access during procedures. 41.3 inches



## Features

Siderails have integrated controls for both patient and caregiver that easily rotate down for patient transfer or access.



### SIMPLE CONTROL

Integrated control panel with clear graphics allow for simple operation by patient and caregiver.

#### Controls:

- 1. Height and Tilt Adjustment
- 2. Trendelenburg Position
- 3. Under Bed Light
- 4. Nurse Call



### SIMPLIFIED CLEANING

High hygienic standard and compact construction and unique features of the bed base reduce cleaning time.

- 1. Seamless upholstery.
- 2. Simple mattress removal for easier cleaning access.
- 3. Epoxy powder coating paint for better and corrosion resistance.

The smooth design angles, column design and open architecture improve safety by reducing pinch points and sharp edges.





Trendelenburg



CPR

#### **FEATURES**

- 1. Trendelenburg
- 2. CPR
- 3. Battery backup
- 4. Easy to clean and disinfect
- 5. Innovative column design
- 6. Easy to access central brake pedal
- 7. Steering castor and handles for better maneuverability
- 8. Removable headboard for anesthesiologist
- 9. Nurse Call
- 10. Tiliting holders for IV pole

# Positioning options for



Lithotomy



Side-lying



Semirecumbent position with the feet rested on the leg rests



Upright



Semirecumbent position with the feet rested on the foot section



Labor bar

# every patient



All fours position



Half sitting using leg holders



Half sitting using foot section



Epidural procedure



Perineal repair



Half sitting with partner assistance

## Accessories







Stainless steel bowl with castors



Telescopic infusion stand



Upholstered support bar



Mattress overlays



Doctor's chair



Doctor's chair, ergonomic

(Ask your local LINET Account Executive about available Doctor Chair colors.)

### STANDARD EQUIPMENT

- Standard leg rests with leg straps (pair)
- Battery backup
- Sliding foot section including the upholstery
- Seamless mattress with bacteriostatic leather
- Electric Trendelenburg and manual CPR
- Infusion stand holder on both sides of the bed
- Siderails with controls
- Fixed panel with controls
- Sliding stainless fluid basin, 4.5 liter

- Hand grips (pair)
- Double castors, 6 in
- Central brake
- Equipotential coupling
- Configurable nurse call cable
- Under bed light
- Pillow
- Removeable headboard
- Front holder

#### TECHNICAL PARAMETERS

Total length of the lying area	83.5 in
Minimum length	61.8 in
Full width	38.6 in
Castor size	6 in
Minimum height	23.6 in
Maximum height	41.3 in
Back section	-10°/+70°
Seat section	0°/+18.5°
Trendelenburg position	0°/-10°
Leg rests, vertical	0°/+135°
Leg rests, horizontal	0°/+60°
Foot section adjustment angle	0°/-22°
Maximum patient load	529 lb (240 kg)
Overall product weight	540 lb (245 kg)



## AVE 2

- Freedom to change a birthing position
- Yes to natural childbirth
- Maximum safety for baby and woman in labor



- Intuitive and easy operation
- Simplified cleaning
- Minimum physical effort



šta∏

anagement

- Exceptional modern design
- Increased level of care
- LDR and LDRP concept



