


Protevo APC

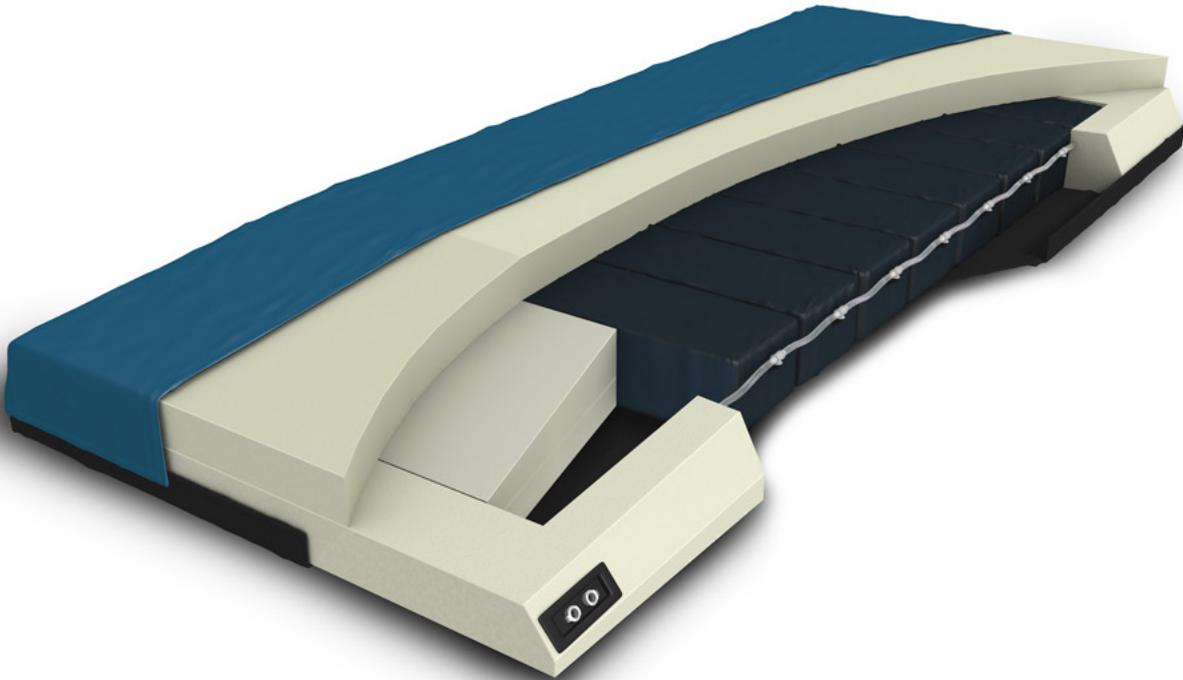


CONVERTIBLE PREVENTION AND TREATMENT

Protevo APC is an advanced pressure redistribution mattress that provides flexible therapy to help meet the wound care needs of your patient. Used without the control unit, Protevo APC provides pressure redistribution through self-adjusting air cells. Easily convert to an active surface with optional control unit to provide alternating pressure therapy.

- 8 air cells absorb and displace patient weight.
- Provides static pressure redistribution without a pump.
- Alternating pressure therapy provided in 10 minute cycles when connected to a pump.
- Dedicated heel zone with 5° slope.
- Heel zone features 2 air cylinders and viscoelastic foam pillow.
- Comfort layer spans entire length and width of mattress.
- Non-slip mattress bottom for patient safety.
- Multi-directional stretch top cover with full zipper for easy replacement.
- Side bolster foam for support.

Protevo APC at a glance

The Protevo APC mattress can be used with any flat deck frame. When used with a LINET® bed, the Protevo APC mattress works in combination with the ErgoFrame® design to help maintain patient position, neutralize shear and minimize pressure in the sacral area.

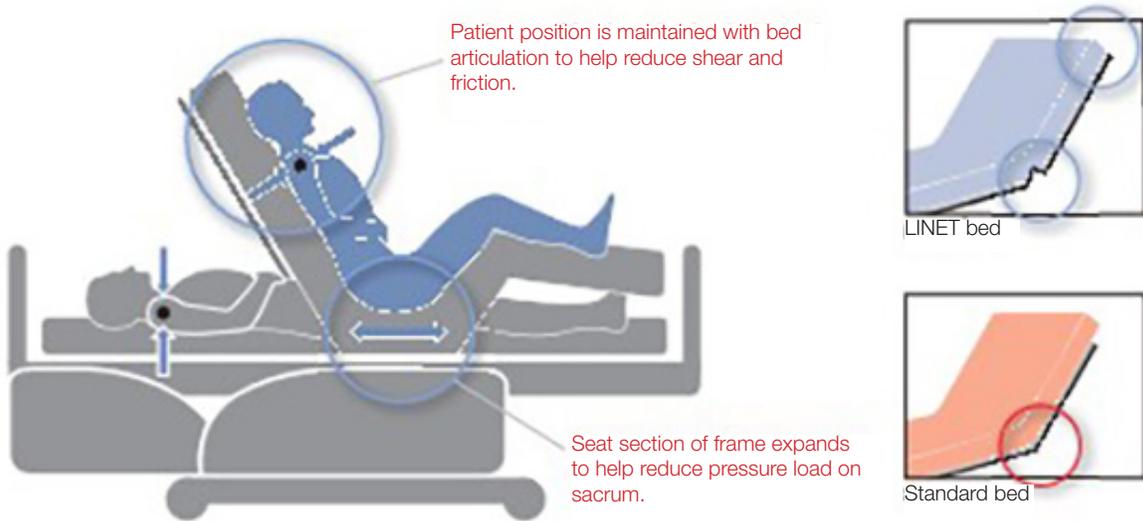

TECHNICAL PARAMETERS*

Width	35 - 35.5 in	Weight Capacity	500 lb
Length	80 - 84 in	Air Foam Cells	8
Height	7 in	Flammability	Meets Cal TB 129, 16 CFR 1633
Heel Slope	5°		

*Width, length and height depends on bed frame model. Compatible LINET® bed frames include Eleganza Series and Multicare® Series.

LINET Americas, Inc.

Phone: 877-815-9895 | E-mail: info@linetamericas.com | www.linetamericas.com

© 2021 LINET Americas, Inc. | 060121
 LINET®, ErgoFrame®, and Multicare® are registered trademarks of LINET Group SE