

# Sentida



## STANDARD FEATURES

- SafeFree Side Guards: 4-stage adjustment
- Adjustable Bed Length: 3 positions
- BiNetic Patient Surface: double retraction
- 3-Stop Strategy: Fall Prevention (sleep), Mobilization (ingress/egress) and Nursing (injury reduction)
- Automatic Positions: Comfort Chair and Semi-Fowler
- Braking at foot of bed
- Warranty: 3 year parts and labor

## SURFACE OPTIONS

Choose from three different patient surface options:

- Aero: Wire mesh surface ensures good aeration.
- Easyclean: Surface with synthetic mats for easy cleaning.
- Comfort: Promotes bedding comfort and makes a valuable contribution to pressure injury prevention with 55 spring elements.


## SPECIFICATIONS

Safe working load	500 lb
Frame height - low	10 in.
Frame height - high	31.5 in.
Frame width (fixed)	35 in. or 39 in.
Frame length (adjustable)	78.75 in., 82.75 in., 86.75 in.
Double Retraction	Yes
Semi-Fowler	30 deg.
Comfort Chair	Yes
Color Options	14 colors
Design Option (Head/Foot)	7 styles



## Easy, Ergonomic

The automatic positions of Sentida, adjustable at the push of a button, increase comfort when lying in bed and support pressure injury prevention for the benefit of the resident while reducing the strain on nursing staff. The automatic comfortseating brings the bed into a pre-selected seating position for watching TV and reading comfortably.

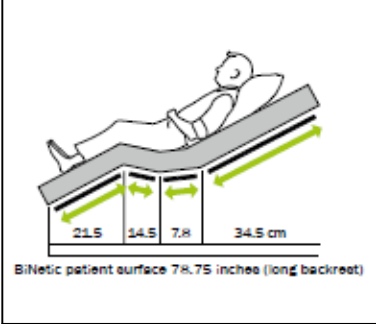

## Effective pressure injury reduction

The ergonomic division of the patient surface helps reduce pressure and therefore contributes to pressure injury reduction. The dimensions of the 4-part patient surface of our nursing home beds are based on the average size and height of a person. With the double retraction, not only the backrest, but also the leg rest are located further apart from the seat section so that the surface is enlarged.




## Flexible and time-saving

All standard wissner-bosserhoff nursing home beds are equipped with an integrated bed extension. With the express bed extension, Sentida - with divided side guards - can be adjusted flexibly and without tools to three different bed lengths. The requirements of the standard IEC 60601-2-52 are complied within each configuration.



## Watching television, reading or eating breakfast in bed with automatic comfort seating

The automatic comfort seating brings the bed into a preselected seating position for watching TV and reading comfortably, for eating or communicating at eye level. Mobilization using the lifter is effectively supported in this position. In addition, this position facilitates breathing and digestion. It trains spatial orientation and stimulates circulation for residents confined to bed. The double retraction provides additional relief of pressure in the sacral area.

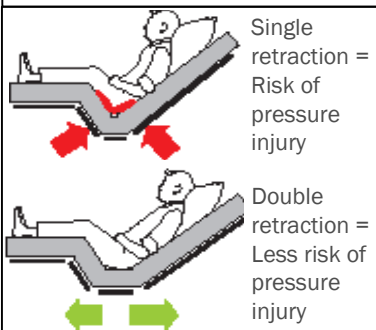

## Effectiveness of double retraction confirmed

For years now, wissner-bosserhoff has been pursuing a practicable concept for reducing the risk of pressure injuries: The double retraction as a component of the BiNetic patient surface. With the double retraction, not only the backrest, but also the leg rest is located further apart from the seat section so that the surface is enlarged. Pressure measurement results show that there is less pressure with double retraction than with single retraction and significantly less pressure than with the standard patient surface.