

Eleganza 3



High-functioning, acute care bed ideal for general medical surgical settings.

Utilize the easy-to-operate and intelligent design to increase caregiver efficiency and patient comfort.



Ergoframe®

Unique surface platform extends space in the sacral area when positioning the bed to reduce pressure.



Early Mobilization

Incorporate standard features for active patient participation during mobilization.

Eleganza 3 at a glance

1 Ergoframe®
Reduces shear and friction forces caused by patient migration to the foot of the bed

2 Mobi-Lift®
Built-in sit-to-stand device for patient to raise and lower the bed

3 Integrated LED Light
Bright LED light in patient handset

4 i-Drive
Automatic 5th wheel enhances maneuverability during transport

5 Easy Access
Removable flat deck for easy cleaning and maintenance


6 Custom Colors
Choose from a variety of footboard and siderail colors

7 One Touch Controls

- Cardiac Chair
- CPR
- Autocontour
- Examination position
- Bed egress position
- Vascular knee adjustment
- Trendelenburg
- Lock-out functionality

8 Hands-Free Foot Controls
Available on both sides of the bed:

- Bed height
- Exam position
- Left and right turn assist
- Brakes

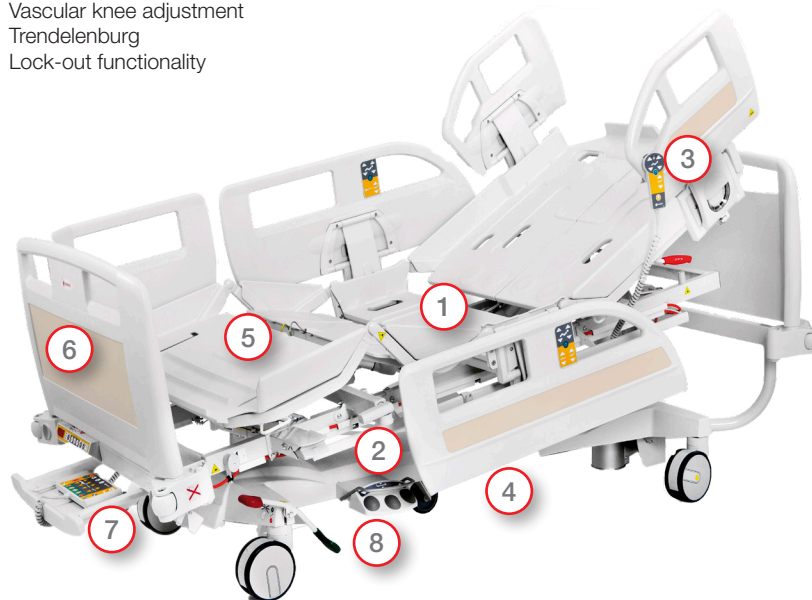

TECHNICAL PARAMETERS

Safe working load	550 lb	Frame weight	315 lb
Frame height - low	15.5 in	Frame width	41 in
Frame height - high	30.5 in	Trendelenburg	± 16°
Frame length extended	97.7 in	Max knee angle	35°
Frame length retracted	84.4 in	Max HOB angle	70°

wissner-
bosserhoff



Members of LINET Group

LINET Americas, Inc.

Address: 10420-R Harris Oaks Blvd., Charlotte, NC 28269

Phone: 877-815-9895 | Fax: 704-248-5655 | E-mail: info@linetamericas.com | www.linetamericas.com

© 2020 LINET Americas, Inc. | 040720

LINET®, Ergoframe®, and Mobi-Lift® are registered trademarks of LINET Group SE