

Multicare[®]

High acuity critical care bed



STANDARD FEATURES

- Platform-based, programmable Automatic Lateral Therapy[®]
- Lateral tilt turn assist
- Mobi-Lift[®] sit-to-stand
- i-Brake[®] automatic locking brakes
- i-Drive[®] automatic 5th wheel
- Dual touch screen displays with built in user's guide
- Scale with weight history
- Multi-zone bed exit
- Hands-free foot controls
- One touch cardiac chair
- One touch 30 degree HOB
- One touch egress position
- Electric bed extension
- Patient handset with flashlight
- Quiet DC motors with battery back-up
- Accessory outlets
- Dual 6" Tente[®] casters
- Accessory outlets standard

SURFACE OPTIONS

Multicare is designed with open architecture to allow flexibility to choose a surface that meets your specific wound care needs. Regardless of the surface, our ErgoFrame[®] positioning system moves in harmony with the surface to help reduce pressure in the sacral area as the bed articulates.

- OptiCare integrated air mattress with constant low pressure and microclimate management
- Protevo[™] advanced pressure redistribution surfaces with convertible options
- Ergozone[™] prevention mattresses

SPECIFICATIONS

Safe working load	550 lbs.
Frame height low	17.3 in.
Frame height high	32.4 in.
Frame width	41.3 in.
Frame length extended	93.3 in.
Frame length retracted	84.6 in.
Frame weight	496 lbs.
Lateral platform tilt	±30 deg.
Trendelenburg	13 deg.
Reverse Trendelenburg	16 deg.
Max knee angle	35 deg.
Max HOB angle	70 deg.

*A lower egress height may be obtained by utilizing the tilt function.




Automatic Lateral Therapy®

- Platform-based rotational therapy
- Simple programming of turn angles and pause times
- Rotate patients up to 30 degrees
- Positioning pads available




Mobi-Lift® Sit-to-Stand

- Built-in sit-to-stand device
- Handles on both sides of the bed
- Simple controls allow raising and lowering of bed
- Allows active participation from patient during mobilization

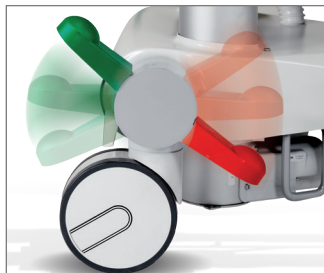

i-Brake®

- Intelligent brakes automatically engage after 60 seconds when bed is plugged into power
- Brake and steer pedals on all 4 wheels

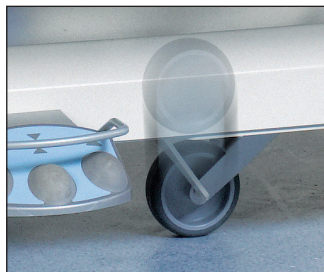

i-Drive®

- Automatic 5th wheel
- Enhances maneuverability of bed during transport
- Retracts when not in use



C-Arm Compatible

- Column construction provides open under bed clearance to accommodate a majority of C-arms
- Radiolucent deck



Multiboard®

- LCD touch screens on both sides of the bed with icon driven controls
- Easy to use interface
- Built-in user's guide
- Advanced patient weight system
- Digital HOB angle indicator
- Digital trendelenburg angle indicator




Lateral Tilt Turn Assist

- Unique platform-based turn up to 30 degrees
- Hands-free foot controls for ease-of-use (15 deg.)
- Open architecture mattress platform




X-Ray Cassette Holder

- Radiolucent deck
- Simplified chest x-ray without disturbing patient

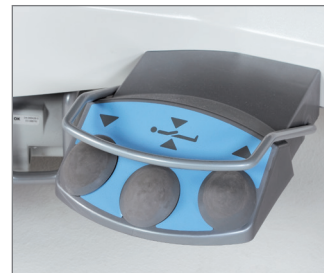

Hands-Free Foot Controls

- Hands-free control of bed
- Bed height
- Exam position
- Left and right turn assist
- Both sides of the bed

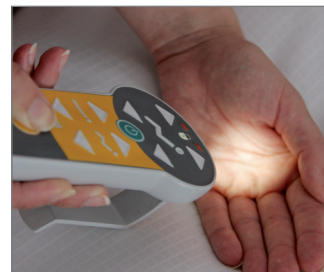

Integrated LED light

- Bright LED light for night time exams
- Integrated into patient handset



Ergoframe®

- Helps reduce patient slide to foot of bed
- Auto-contour button raises head and knee simultaneously to maintain position and increase comfort

One Touch Controls

- 30 degree HOB
- Cardiac Chair
- CPR
- Auto-contour
- Examination position
- Bed egress position
- Vascular knee adjustment
- Electric bed extension