


Protevo GTE

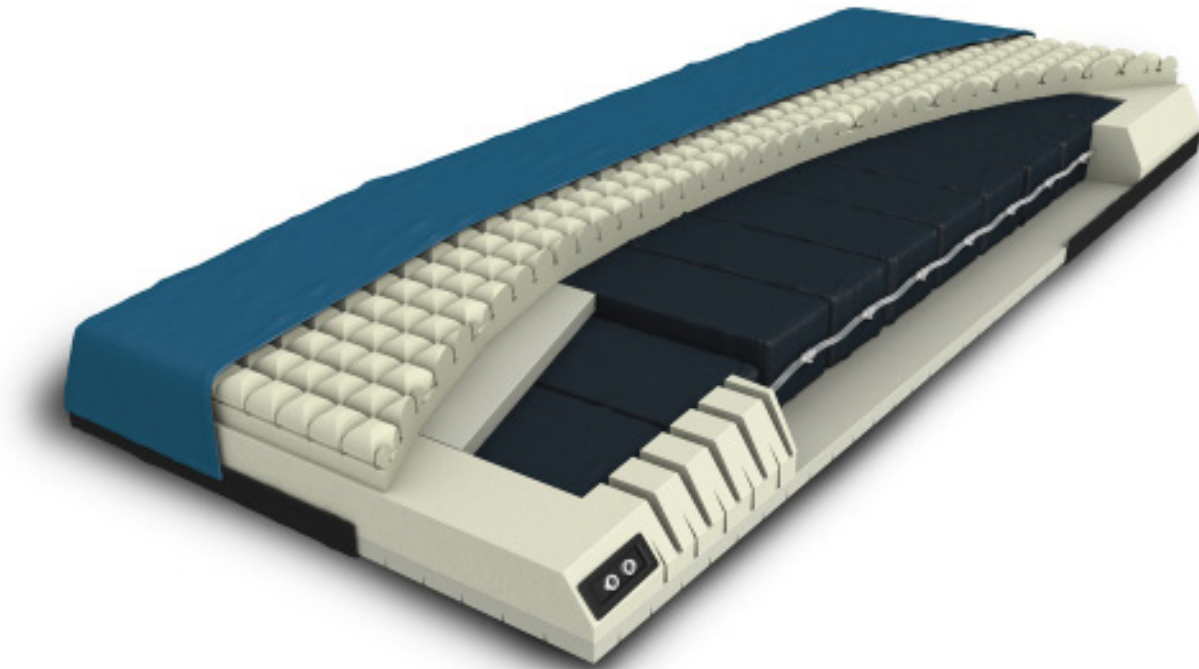


CONVERTIBLE PREVENTION AND TREATMENT

Protevo GTE is an advanced pressure redistribution mattress that provides flexible therapy to help meet the wound care needs of your patient. Used without the control unit, Protevo GTE provides reactive pressure redistribution through self-adjusting air cells. Easily convert to an active surface with optional control unit to provide alternating pressure therapy.

- Increase patient comfort and assist in the reduction of pressure injury incidents with the 3-dimensional, sculpted envelopment layer.
- 10 air cells absorb and displace patient weight.
- Provides advanced reactive pressure redistribution.
- Alternating pressure therapy provided in 10 minute cycles.
- Dedicated heel zone with 5° slope.
- Heel zone features 2 air cylinders and viscoelastic foam pillow.
- Comfort layer spans entire length and width of mattress.
- Non-slip mattress bottom for patient safety.
- Multi-directional stretch top cover with full zipper for easy replacement.
- Side bolster foam for support.

Protevo GTE at a glance

When used with a LINET® bed, the Protevo GTE mattress works in combination with the Ergoframe® design to help maintain patient position, neutralize shear and minimize pressure in the sacral area.

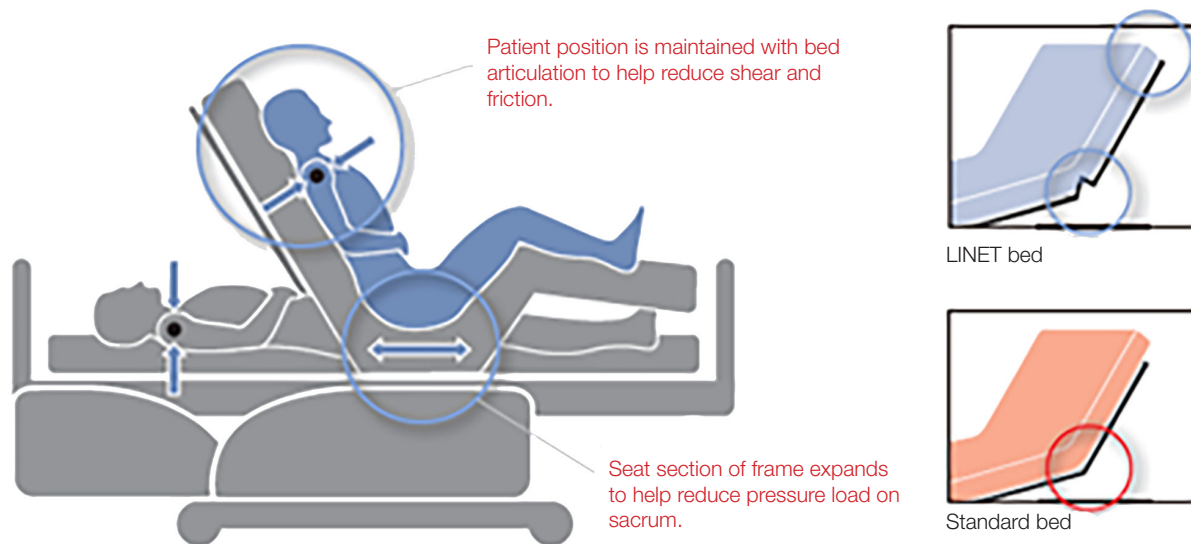

TECHNICAL PARAMETERS*

Width	35 - 35.5 in	Weight Capacity	600 lb
Length	84 - 85.5 in	Air Foam Cells	10
Height	7 - 8 in	Flammability	Meets Cal TB 129, 16 CFR 1633
Heel Slope	5°		

*Width, length and height depends on bed frame model. Compatible LINET® bed frames include Eleganza 3, Eleganza 5, and Multicare® and Multicare LE



Members of LINET Group

LINET Americas, Inc.

Address: 10420-R Harris Oaks Blvd., Charlotte, NC 28269

Phone: 877-815-9895 | Fax: 704-248-5655 | E-mail: info@linetamericas.com | www.linetamericas.com

© 2020 LINET Americas, Inc. | 040820
LINET®, Multicare® and Ergoframe® are registered trademarks of LINET Group SE