


Protevo ST

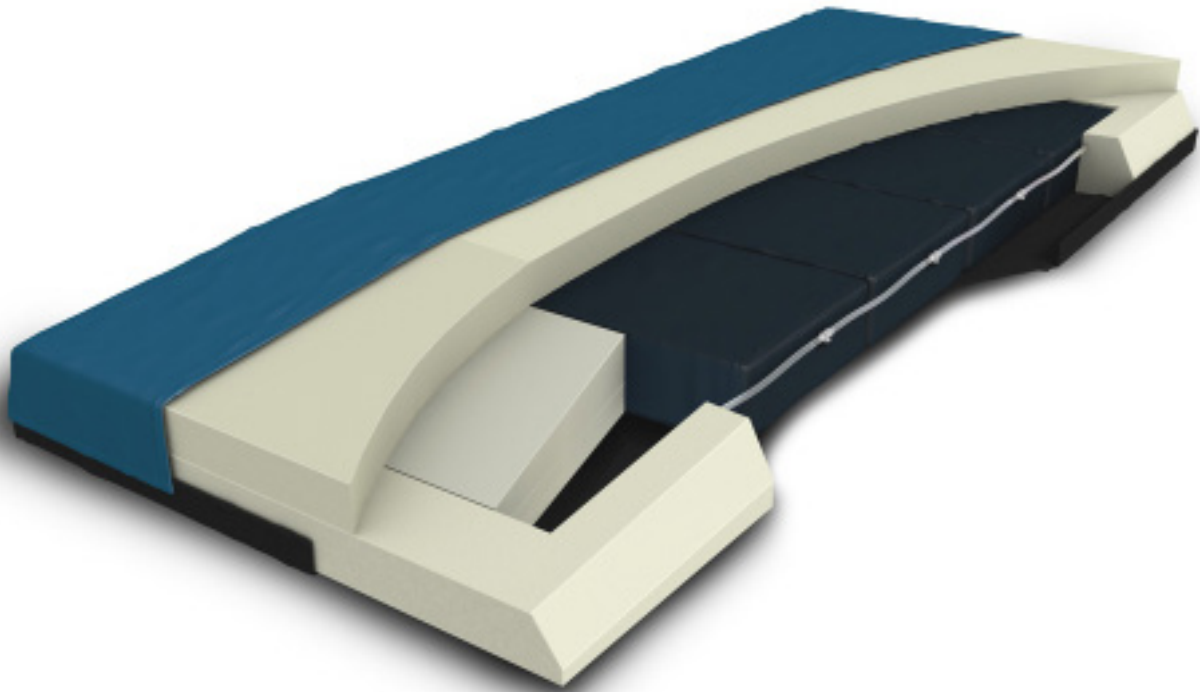


COMFORTABLE, RESTFUL RECOVERY

Protevo ST is an advanced pressure redistribution mattress offering optimal envelopment through self-adjusting air cells. Simple, economical wound care designed for the comfort of your patient.

- 4 air cells beneath comfort layer to displace patient weight.
- Self-adjusting, advanced pressure redistribution without a control unit.
- Dedicated heel zone with 5° slope.
- Non-slip mattress bottom for patient safety.
- Side bolster foam for support.
- Comfort layer spans entire length and width of mattress.

Protevo ST at a glance

The Protevo ST mattress can be used with any flat deck frame. When used with a LINET® bed, the Protevo ST mattress works in combination with the Ergoframe® design to help maintain patient position, neutralize shear and minimize pressure in the sacral area.

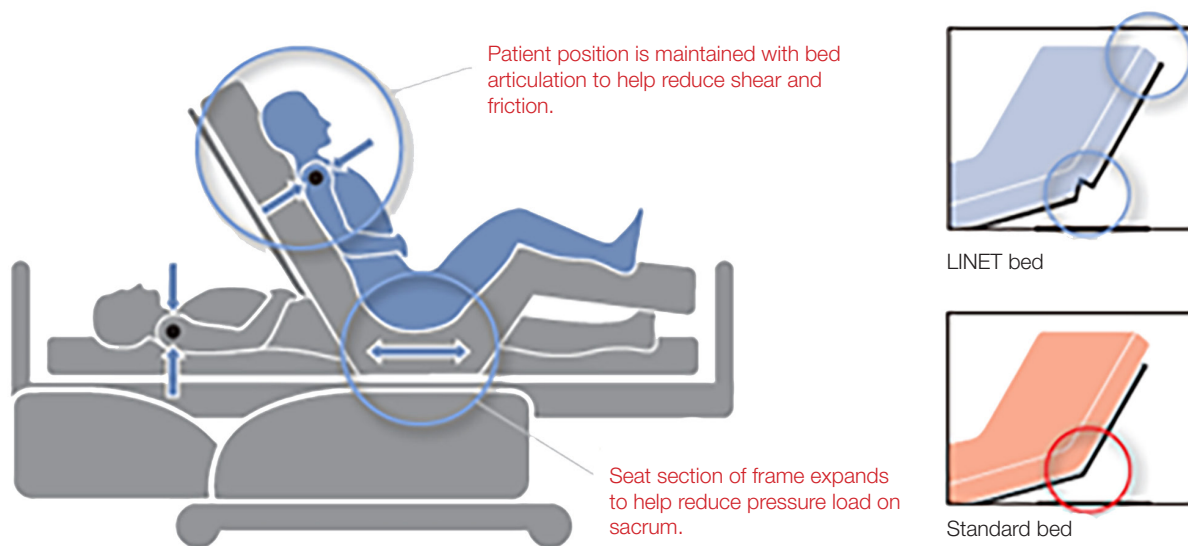

TECHNICAL PARAMETERS*

Width	35 in	Weight Capacity	500 lb
Length	84 in	Foam Cells	1
Height	7 in	Air Cells	4
Heel Slope	5°	Flammability	Meets Cal TB 129, 16 CFR 1633

*Width, length and height may depend on bed frame model. Compatible LINET® bed frame includes Eleganza 3.



Members of LINET Group

LINET Americas, Inc.

Address: 10420-R Harris Oaks Blvd., Charlotte, NC 28269

Phone: 877-815-9895 | Fax: 704-248-5655 | E-mail: info@linetamericas.com | www.linetamericas.com

© 2020 LINET Americas, Inc. | 040820

LINET® and Ergoframe® are registered trademarks of LINET Group SE